

Updating EASYSense Vu+ Firmware


Occasionally Data Harvest may release updated firmware which will contain improvements or new features. Updates will be made available from the [Data Harvest website](#).

To update the firmware of your **EASYSense Vu+** logger:

Step 1. Copy the Firmware Update for Vu+ onto the hard disc of your computer and run the .exe.



Step 2. Connect Vu+ to the computer using the USB cable.

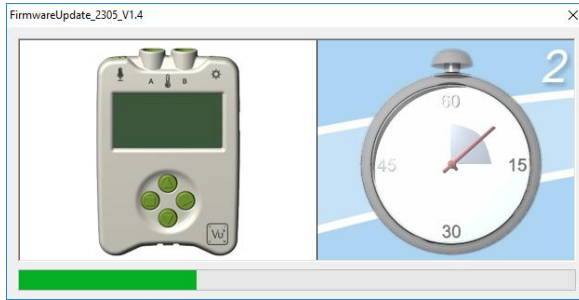


Step 3. Wait for a few seconds.

A progress bar will indicate Vu+ is being programmed.



A **tick** will indicate if the process was successful. Disconnect Vu+.

If you wish to repeat this process, connect another Vu+ and click on the update window.

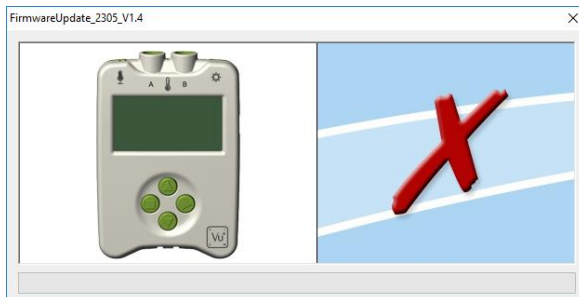

A **cross** will indicate that an error occurred. Click on the update window to go back to the start screen and re-program.

Step 4. Exit the program.